



WELL Communications


Vico Care

by Haelvoet



HAEVOET NV
L. BEKAERTSTRAAT 8, 8770 INGELMUNSTER, BELGIUM
T +32 (0) 51 48 66 95 F +32 (0) 51 48 73 19
INFO@HAEVOET.BE WWW.HAEVOET.BE

HOSPITAL - NURSING HOME FURNITURE

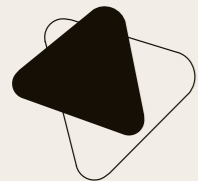

Vico

technology you can rely on



001	Mattress platform • For a comfortable lying position
	You can choose between 2 different mattress platforms. The Vico is standard equipped with well-ventilating metal mattress platform profiles. The optional solid laminate mattress platform sections (*) can be easily removed and cleaned.
002	Mattress holders • Purposefully present
	The mattress holders prevent the mattress from being moved towards the bed sides. They can be also used as a handle to adjust the backrest and knee rest manually.
003	Sideguards (*) • Safety has no price
	The Atmosphere sideguards protect the patient over 2/3 of the bed length, whereas the Pollux and Trix sideguards offer protection over the entire length.
004	Cases • Optimal care
	Two cases at the head end can be used to install an IV rod, a lifting pole or other accessories. Thanks to this, the bed can be equipped in accordance with the user's specific needs.
005	Wooden surrounds of the mattress platform frame (*) • An extra trump
	The wooden surrounds on the mattress platform frame give the bed a homely look.
006	Operating lever • Central and within reach
	This lever activates a gas spring, so that the desired mattress platform section can be adjusted.
007	Cover cap (*) • Hygienic
	The removable rounded synthetic cover cap makes it easy to clean the carriage.
008	High/low system • Highly stable
	The high/low supports guarantee excellent bed stability at all heights. The bed can be easily adjusted between 41 cm and 80 cm (with castors \varnothing 125 mm), which makes it easier to get in and out of the bed or to treat a patient.
009	High/low pedal • Easy to operate
	A hydraulic pump guarantees the high/low movement. It can be adjusted by means of a pedal.
010	Central brake (*) • Activated with one single movement
	The castors can be centrally braked by means of a pedal within easy reach. This simplifies the use of the bed by the nursing staff.
011	Castors • Superb maneuverability
	Our extensive gamut of high-quality design castors guarantees a highly mobile and maneuverable bed.
012	Bed panels • A harmonious duo
	The head and foot panels are available in 4 different versions with stylish wood designs.

Haelvoet nv
Edition 2010
Haelvoet reserves the right to modify the product.
Small colour deviations are possible.


* option for Vico Care

Comfort for everyone

The easy operation and efficient user-friendliness of the Vico Care bed support the nursing staff during their daily work.


- Activations**

The operating lever is user-friendly and easy to activate. The high/low pedal is within easy reach on both sides of the bed.

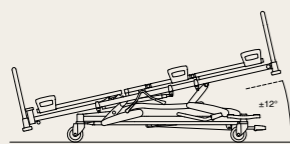


- (Reverse) Trendelenburg position**

The (Reverse) Trendelenburg position (*) can be activated by means of the red pedal at the foot end. When the pedal is stepped on, the bed inclines fluently until the desired position is reached.

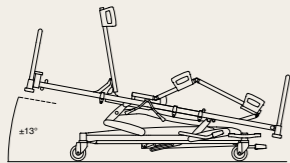

- Seating or relax position**

The 3-section or 4-section (*) mattress platform offers ample adjustment possibilities and can be easily activated on both sides of the bed. Thanks to this, the desired seating position – or relax position in case of a 4-section mattress platform (*) – can be easily obtained.

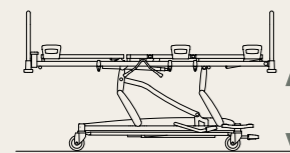

- High/low position**

The high/low pedal enables you to adjust the bed height on both sides of the bed, without needing your hands. The highest position of the bed is 80 cm (castors \varnothing 125 mm). Daily care and treatments can be performed without having to strain the nurse's back.



- Retracting backrest**

A retracting backrest increases the resident's lying comfort considerably. While inclining the backrest upwards, it is automatically moved backwards as well. By doing so, pressure on the pelvis is relieved, whereas the bedsore risk is reduced substantially.




Mobile and flexible


The daily use is made more comfortable, thanks to an easily movable and manoeuvrable bed.

- Central brake**

Both sides of the foot end are equipped with a design brake pedal (*). This pedal activates the 4 castors simultaneously, and can be put in 3 different positions.

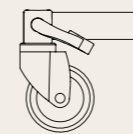


- 4 braked castors**

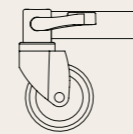

- 4 unbraked castors**

This position makes it possible to move the bed in all directions.

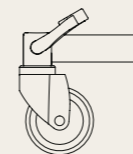

- 1 fixed directional castor**

This castor makes it easier to move the bed in a straight line.



User-friendly

- Mechanical**


- A minimum of maintenance**

A maintenance-free hydraulic pump and maintenance-free gas springs make this bed a durable investment.

- Removable cover cap**

- Hygienic**

Thanks to the rounded synthetic cover cap (*), the surfaces that are difficult to reach can be cleaned easily.

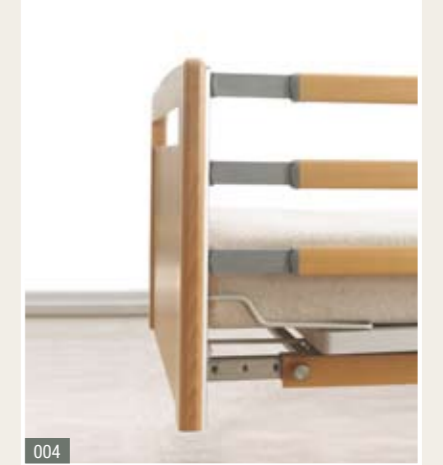

- Bed extension**

A nursing bed has to fulfil everybody's wishes, irrespective of the resident's length. The integrated bed extension makes it possible to extend the bed with 20 cm.

Two different methods are possible to extend the bed. You can choose between a telescopic bed extension combined with the extendible Trix sideguards or a permanent bed extension (of 10 or 20 cm) when using the long Pollux sideguards.

When combined with the optional foot extension (*), a complete mattress platform functionality is guaranteed.

- No current**

The Vico Care bed is completely mechanical. Thanks to this, the bed can be used everywhere without any problems.

001 The operating lever.

002 The high/low pedal is installed on both sides of the bed, and is within easy reach of the nursing staff.

003 The (Reverse) Trendelenburg pedal (*).

004 When the bed has been extended, the telescopic Trix sideguards (*) offer complete protection over the entire bed length.


005 The single brake pedal is standard, in case the central brake mechanism (*) is chosen. The double brake pedal and the directional castor are optional.

006 The removable cover cap (*) makes it easy to clean the bed.

* option for Vico Care

Safety galore


Our different sideguards offer complete resident safety in each bed position. Each type has its specific advantages for both the user and resident. The choice is completely yours.

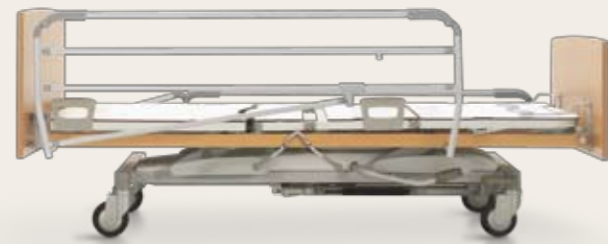

The Atmosphere sideguards (*) are one unit and can be lowered by means of an easy locking mechanism.



The Pollux sideguards (*) consist of 2 strong wooden siderails. A guaranteed distance between the two siderails makes sure that the Pollux sideguards comply with the strictest norms. The sideguards can only be lowered by using your both hands. This safety system has been developed to prevent the sideguards from being lowered unintentionally.

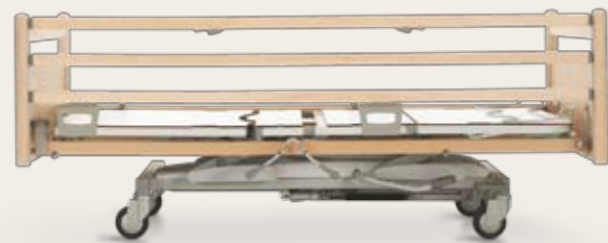

The Trix sideguards (*) consist of 3 siderails, offering maximum safety over the entire bed length. An integrated telescopic system makes sure that the Trix sideguards – even when the bed is extended – cover the bed sides over the entire length.

007 When thicker anti-bedsore mattresses are used, an optional 4th siderail () can offer 13 cm of extra protection. This protection is completely integrated in the Trix sideguards, so no additional accessories are needed.*

008 The sideguards can only be lowered by using your both hands. This safety system makes sure that the sideguards cannot be lowered unintentionally.

* option for Vico Care












007



008

Accessories

The accessories we offer can be harmoniously connected with the existing cases and accessory holders. Haelvoet makes sure that the accessories are at least as good as our nursing beds, so that you can always have absolute certainty that your product is top quality.

001	Lifting pole		The lifting pole is made of a bent thick-walled steel tube. It is equipped with a handle and a ribbon that can be manually adjusted. An automatic roll-up system is optionally available as well.	006	IV rod		The IV rod consists of a fixed and a height-adjustable part. The latter is equipped with four hooks. The IV rod is also available in a non-adjustable version with 2 hooks.
002	IV holder for lifting pole		The IV holder for the lifting pole has 4 hooks and can be easily installed on the lifting pole.	007	Raising aid		The ergonomic raising aid offers the resident optimal support when getting in and out of the bed.
003	Dinner tray for sideguards		The wooden dinner tray can be installed on all sideguard types. It can be easily hung onto the foot panel.	008	Urine bottle holder		The urine bottle holder is installed in the 2 synthetic accessory hooks that are optional for the Vico bed. You can choose between 2 models.
004	Accessory bracket		The accessory bracket can be used in case of specific customer needs. The bracket is equipped with 4 hooks and has a length of 310 mm. It is put in the optional accessory hooks.	009	Cover for sideguards		The removable cover is made of imitation leather and can be easily put around the sideguards. This cover is especially recommended to protect restless residents.
005	Bowl holder		Two bowls can be placed into the bowl holder. The holder is put in the optional accessory hooks.				

Elegant design

Everybody cares about the health of residents. Personal independence and design play a key role in this concept. That's why the most modern techniques and a progressive functionality have been integrated in the Vico Care bed design.

- Bed panels**

Our bed structure makes it possible to choose between 4 different panel types. Each panel is characterised by its specific lines and sharp eye for detail. Thanks to the absence of visible fastening bolts, the panels don't only have a stylish look, but they are also very easy to maintain and clean.



009 The matching colours, homelike materials and functional design improve the resident's health.



Type 1: The standard bed panel **Midea** has a modern look.



Type 2: The timeless bed panel **Hypnos (*)** can be easily integrated in various interiors, thanks to its round solid posts.



Type 3: The creative bed panel **Felina (*)** is characterised by a slightly bent handle.


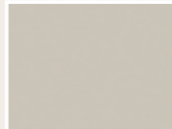






Type 4: The bed panel **Valentino (*)** is characterised by rigid lines.

- Finishing and quality**

A whole range of colours from our standard collection helps you to integrate this bed perfectly in your institution.

Melamine

 snow white code 8685BS	 cashmere code 5981BS	 classic beech code 0381BS	 light select walnut code K008PW
 birch code 1715BS	 blonde liberty elm code K017PW	 wengé code 0854BS	 louisiana wengé code 9763BS

Solid wood

 nature varnish code 25	 white varnish * code 121	 washed oak * code 170	 havana * code 056
 puce dark oak * code 020	 wengé * code 125	 cashmere* code 149	 white mat lacquer * code 300

Samples from our stained colours are available on request.

Other solutions for your interior decoration are possible on request.




This PEFC certificate guarantees that Haelvoet uses wood of sustainably managed forests. A sustainable forest management is environmentally and economically viable, as well as socially beneficial. By buying our PEFC products, you contribute to the future of our environment.

Wood is a natural product. Small colour deviations are possible.



Vico Care Summary

• Dimensions and performances (standard bed)

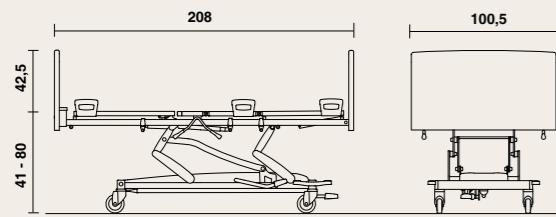

- width: 100,5 cm, length: 208 cm;
- height variation of the mattress platform:
 - ± 41 cm - ± 80 cm (with castors Ø 125 mm),
 - ± 40 cm - ± 79 cm (with hidden castors);
- clearance height: 100 à 150 mm (depending on the height adjustment of the mattress platform);
- safe working load: 200 kg;
- weight (without accessories): 115 kg;
- mattress dimensions: 195 x 85 cm, thickness: 12 to 15 cm.



Haelvoet makes sure that its excellent price-quality ratio is not at the expense of the environment. Our company considers durability – in all possible meanings – of paramount importance. All our products have a long lifespan and are delivered with a 10-year-guarantee on construction errors.



Choosing Haelvoet is choosing an honest product, thanks to our environment-friendly production processes, a maximum use of recyclable materials and our PEFC certificate.

CE The Vico Care nursing bed is made in accordance with ISO 9001:2000 and complies with the European Medical Devices Directive 93/42/EEC and all other applicable European harmonised norms:

EN 60601-1 Medical electrical equipment - Part 1: General requirements for safety, except for the electrical regulations.

EN 60601-2-52 Medical electrical equipment - Part 2-52: Particular requirements for the basic safety and essential performance of medical beds, except for the electrical regulations.

EN ISO 14971 Application of risk management to medical devices.

	Standard Vico		Options Vico			
MATTRESS PLATFORM						
	• 3-section steel mattress platform	• Telescopic bed extension up to 20 cm	• 4-section steel mattress platform	• 3-section mattress platform with removable solid laminate sections	• 4-section mattress platform with removable solid laminate sections	• Foot extension up to 20 cm
CASTORS						
	• Levina castors Ø 125 mm, with separate brakes		• Cover design castors Ø 125 mm, central brake, with or without directional castor	• Cover design castors Ø 150 mm, central brake, with or without directional castor	• Integral castors Ø 150 mm, central brake, with or without directional castor	• Hidden castors Ø 50 mm, with central brake
PANEL TYPE						
	• Midea panel (only with Atmosphere and Pollux sideguards)		• Hypnos panel	• Felina panel	• Valentino panel	• Bumpers (not on the Midea panel) • Wall spacer
SIDEGUARDS						
			• Atmosphere sideguards	• Pollux sideguards	• Trix sideguards	• Trix sideguards with 4th siderail
CONTROLS						
	• Operating lever	• High/low pedal	• (Reverse) Trendelenburg pedal			
EXTRAS						
			• Cover cap	• Wooden surrounds of the mattress platform frame	• Decorative strip (front and back profile of the carriage)	• 6 synthetic accessory brackets